




Spectral Correction on PM400

1. Please check that the current firmware (available at https://www.thorlabs.com/software_pages/ViewSoftwarePage.cfm?Code=OPM) is installed on your PM400.
2. Create a .csv file containing your spectral data considering the following points:
 - The spectrum must contain two columns (wavelength, intensity) delimited by a comma.
 - The wavelength values must not contain any decimal places. Therefore, the minimum step size is 1nm, but the intervals can also be larger. The step size does not need to be equal for each step.
 - There must not be any duplicate wavelength entries.
 - The wavelength range must not be larger than that of the sensor. It can be smaller though.
 - The filename must not be longer than 8 characters.
 - The intensity information can be decimal number and the decimal places can differ for each entry.

Please note that creating a .csv file from Excel can lead to a wrong encoding! One option is to just create a txt file containing all data and then rename it to .csv. Otherwise, you can use a text editor to ensure (and edit if needed) that the encoding is set to UTF-8 and not e.g. UTF-8-BOM.

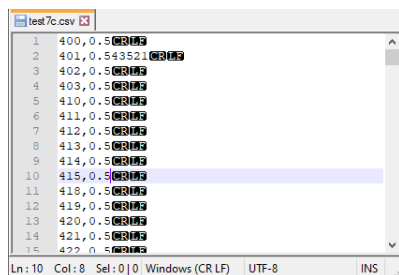

3. Connect your PM400 to your PC via USB and change the Owner to USB. You should see a red icon on the top right:




- ATTEN_CURVES
- LANG
- RECORD
- SPEC_CURVES
- SUPPORT_DOCS

The image shows the screen of a Thorlabs PM4000 Optical Power Meter. The screen displays the 'Spectral Correction' menu. At the top, the Thorlabs logo and model number 'PM4000 Optical Power Meter' are visible. Below this, there are fields for 'Wavelength' (650 nm), 'Power' (0.00000), and 'Power/area' (0.00000). The main menu consists of a grid of options: 605 nm, 632 nm, 633 nm, 660 nm, 785 nm, 830 nm, 850 nm, 1064 nm, TSI01, TEST7C, TEST7g, and 1.3 μm. At the bottom, there is a row of function keys: 'Choice', 'price', 'data', 'Set', 'process', 'test', 'zero', 'cross', 'return', and 'key'. Below the screen, the physical buttons for these functions are visible on the device's faceplate.

The image shows the screen of a Thorlabs PM400 Optical Power Meter. The screen is white with a red border. At the top, the Thorlabs logo and model number 'PM400 Optical Power Meter' are displayed. Below this is a keypad with numbers 1-9, 0, and function keys like 'on', 'off', 'unit', 'mode', 'range', 'zero', 'back', and 'enter'. The screen also shows a date and time display '10/07/2016 10:12:23' and a battery level indicator. At the bottom of the screen, there are four large icons: a wavelength symbol (λ), a triangle, the word 'ZERO', and a list icon.

7. Choose the file you would like to load and click on the folder symbol on the top right.

THORLABS

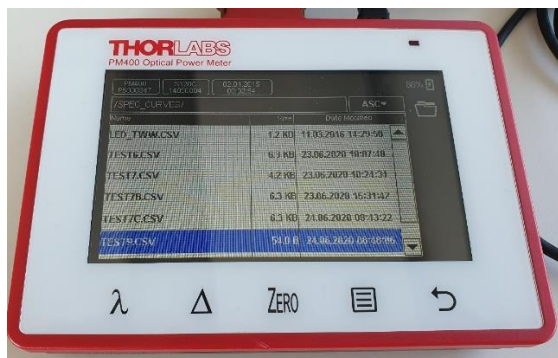

Then click on OK.



8. You should now be able to choose the respective entry in the Spectral Correction menu. PM400 automatically checks, whether the file format and data in the .csv file are correct. If that is not the case, you will get an error notification when trying to choose the spectrum correction entry.

